

The Brilliant Black Reduction Test (BRT)

The BRT Inhibitor Test The BRT MRL-Screening Test

Agar Diffusion Procedure for the Detection of Inhibitors and Medicinal Residues in Milk and Milk Products

Instructions for use


The BRT is carried out as follows:

- Remove the sealing foil from the test plate
- Pipette 0.1ml of the milk sample into one well (a second insertion is also recommended)
- Pipette 0.1ml of the negative control (Inhibitor Free Milk) into one well
- Pipette 0.1ml of the positive control (Penicillin G Standard) into one well
- Seal the test plate with the perforated sealing foil
- Incubate the test plate in a water bath or thermo-block at 65 °C until the negative control has turned yellow (BRT Inhibitor Test, 2 hrs 15 +/- 15 min; BRT MRL-Screening Test, 2 hrs 30 +/- 15 min)
- The milk can then be poured away before evaluation
- The test results are read from the bottom of the test plate

The procedure is the same for test strips or tubes. The test is standardised, and each batch is analysed according to the procedural conditions described above. If the procedure is not followed correctly there can be no guarantee for the accuracy of the results.

The BRT Inhibitor Test is produced according to Commission Decision 91/180/EEC and German Food and Feeding Stuff Law, LFGB, § 64 (before Food and Consumer Goods Law, LMBG, § 35) Method L 01.01-5.

The BRT MRL-Screening Test is produced according to Commission Decision 91/180/EEC and German Food and Feeding Stuff Law, LFGB, § 64 (before Food and Consumer Goods Law, LMBG, § 35) Method L 01.00-11.

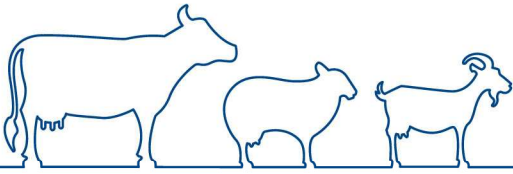

Important notes for execution of the test

Storage and use-by date

The Brilliant Black Reduction Test must be kept cool, and should be stored on delivery at 6-10 °C. Under these circumstances the BRT is usable for at least 3 month (BRT ESL, 6 month; BRT tubes, 9 month). After a period of storage the incubation time may be longer (by about 15-30 minutes, according to length of storage). The detective sensitivity of the product will not however be affected by long storage.

Sample material

The milk samples should be of acceptable quality and should be tested as soon as possible. For short-term preservation a storage temperature of 6-10°C is sufficient. For a longer storage period the material should be frozen at -15 to -30°C to prevent inactivation of any residue. Samples stored in this way should be thawed in a water bath at 45°C, and must be thoroughly mixed before inspection.

Controls/Standards

The negative and the positive controls should always be composed and tested in exactly the same way as the sample material. The influence of the test material on the test system is thereby taken into account.


The **Inhibitor Free Milk** is used in the BRT as a function, or negative, control. The Inhibitor Free Milk is delivered freeze dried and can be pipetted after re-suspending into the test wells. The well containing Inhibitor Free Milk must turn from blue to yellow during incubation.

The **Penicillin G Standard** serves as a positive control. The Penicillin G Standard is delivered freeze dried has to be re-suspended respectively diluted according to the instruction for use and can then be pipetted directly into one of the test wells. The well containing the Penicillin G Standard must remain blue after incubation.

The **Penicillinase solution** can be used for the identification of penicillinase-unstable penicillins. The Penicillinase solution is delivered ready for use and can be used as follows:

Pipette 10µl of the Penicillinase solution into one of the wells. Add 0.1ml of the positive test sample and incubate by the described manner. If the milk sample has been contaminated by penicillinase-unstable penicillins, the well will show a negative reaction after incubation. In the presence of other inhibitors, it will show a positive reaction.

In the case of contamination the solutions or the Inhibitor Free Milk may be rendered unusable or may lead to incorrect results. The controls must therefore always be kept clean and carried out with care.

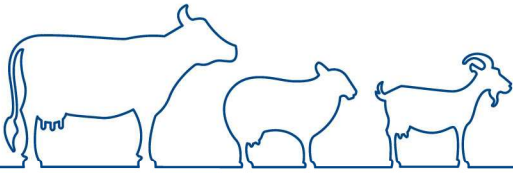

Sealing of test containers

The test plates (strips, tubes) should be sealed with extreme care before being placed in the water bath. If water penetrates the reaction system the test result may be inaccurate.

Incubation

Incubation can take place in a **water bath**, a **thermo-block** or in an **incubation chamber**. The water bath or the thermo-block is fundamentally preferable, in that temperature transfer is more direct and heat distribution more even than in an incubation chamber. In all of these methods it is important that the temperature remain constant and that there be no great discrepancies in the set temperature, either at different times or in different areas.

Incubation temperature

The incubation temperature for the water bath and for the thermo-block is 65 °C. In an incubation chamber the temperature must be set individually according to the type of apparatus, and can be between 65-69°C. The water bath or thermo-block are recommended because of the easier and more reliable regulation of temperature when using these methods. If the incubation temperature deviates from the ideal 65 °C (as measured within the test system), the incubation time will increase considerably (by 30 minutes or more).


Incubation time

The incubation time is dependent on the temperature set as well as on a constantly even distribution of heat. Under ideal circumstances the incubation time is 2 hrs 15 +/- 15 min. The incubation time may, however, increase as a result of factors such as differing storage times. The test must be incubated until the negative control is completely discoloured.

Evaluation

A correct evaluation can only be made if at least one positive and one negative control have been included on the test plate or strip or in one of the test tubes, and both have reacted correctly.

The interpretation of the test-results should be in accordance to the official regulations; these regulations may differ in various countries.

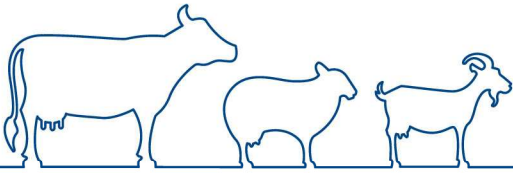

In Germany there are, according to regulations, two possible interpretations of the results:

Inspection tests for **milk quality payment** (inhibitor test in accordance with the official registration of inspection test procedures according to § 64 LFGB (before § 35 LMBG) Method L 01.01-5) are to be evaluated as follows:

- All reaction systems which display at least the colour intensity of the positive control (Penicillin G Standard) are to be evaluated as positive.

Inspection tests carried out in accordance with the German **Foodstuffs Law** (MRL screening test in accordance with the official registration of inspection test procedures according to § 64 LFGB (before § 35 LMBG) Method L 01.00-11) are to be evaluated as follows:

- All reaction systems which display a colour which is clearly to be differentiated from the colour of the negative control are to be evaluated as positive or suspect.

Further information

BRT – Background and Basic Information

BRT – Product Information

BRT – Instructions for Use of Plates

BRT – Instructions for Use of Tubes

BRT – The Correct Incubation

BRT – Use of Preserved Samples

BRT – The Visual Evaluation

BRT – The Photometric Evaluation

BRT – Use of Milk Products

Standards, lyophilised – Product Information

Penicillinase Product Information

Should you have any further questions or require any special information, please do not hesitate to contact us directly.